

Canine Adenovirus Antigen Lateral Flow Assay Kit

Catalog No: E-AD-C065

40T

This manual must be read attentively and completely before using this product.

If you have any problems, please contact our Technical Service Center for help.

Phone: 240-252-7368(USA) Fax: 240-252-7376(USA)

Email: techsupport@elabscience.com

Website: www.vetassay-elab.com

Please kindly provide us the lot number (on the outside of the box) of the kit for more efficient service.

Test principle

This kit applies the principle of Immunochromatography assay. The sample will move together with the colloidal gold marker along the chromatography membrane. If Canine Adenovirus Antigen exist in the samples, it will combine with the colloidal gold marker and the Antibody in the detection line, then it will show a color. Otherwise, it will not show the color reaction.

Kit components

Item	Specification
Detection Card (with a pipette)	40T
Cotton Swab	1 package
Sample Diluent	40 vials
Manual	1 copy

Notes

1. FOR RESEARCH USE ONLY. Please read the manual carefully before use, changes of operation may result in unreliable results.
2. Do not use product out of date or in a broken aluminum foil, it is disposable and cannot be used repeatedly.
3. The detection card should be brought to room temperature before opening after take it out from the refrigerator. The opening detection card should be used as soon as possible.
4. Please do not use but not limited to the following liquids for negative control: water, PBS.
- 5.
6. The tested sample should be fresh and clear. Avoid of using samples of turbidity, polluted, high hemolysis or abnormal viscous.
7. Avoid of touching the chromatography membrane of the sample well and test well.
8. The waste of experiment should be considered as contaminant, and must be properly handled according to the local regulations.
9. Each reagent is optimized for use in the E-AD-C065. Do not substitute reagents from any other manufacturer into the test kit. Do not combine reagents from other E-AD-C065 with different lot numbers.

Storage and expiry date

Storage: Store at 2-30°C. With cool and dry environment, avoid freeze.

Expiry date: expiration date is on the packing box.

Sample preparation

1. **Nasal secretions:** Use cotton swab to collect nasal secretions samples.
2. Insert the swab into the **Sample Diluent** tube and mix the swab until the sample has been dissolved into the assay diluent (Approximately 10 sec).
3. Wait for 1 minute to settle down the large particles.


Assay procedure

Allow all kit components and sample to reach room temperature prior to testing.

1. Tear the aluminum foil bag of the detection card and take out the detection card, and put it on a smooth, clean table.
2. Take the sample with the pipette, add 2~3 drop (about 60 μ L) of sample solution to the sample well vertically and slowly.
3. Incubate for 5 to 10 minutes and then judge the results immediately.

Judgment of result

1. **Negative:** Only the control line region (C) shows a line in the observation well.
2. **Positive:** Both the test line region (T) and the control line region (C) show a line in the observation well.
3. **Invalid:** No line shows in the observation well of the control line region (C).

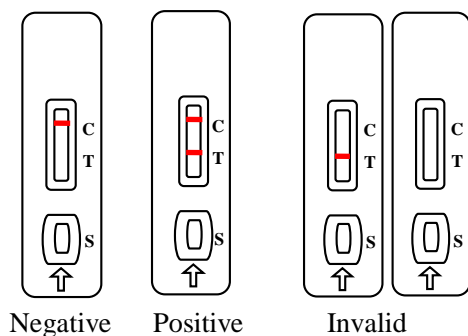

Limitations

1. This kit can only be used for qualitative detection of CAV-Ag in Canine
2. The detection results of this kit are only for reference. For confirmation of the result, please combine the symptoms and other methods of detection, this detection cannot be used as the only criteria for result.